


# SCiZYS

The Future Product for Tissue Analysis

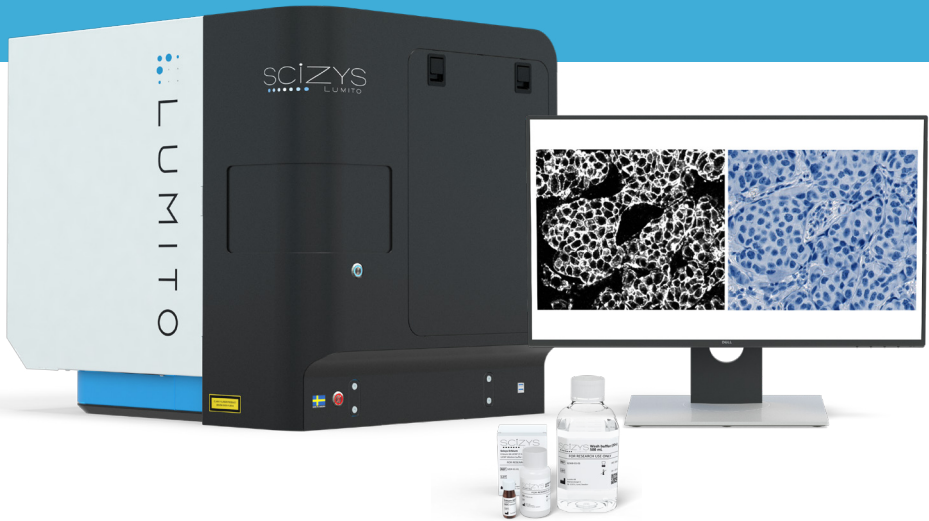

Higher Quality, Higher Accuracy, More Efficient Analysis Processes

More Information  
Without Autofluorescence

Now Available for  
Research Use Only!

LUMITO

# SCIZYS in your Standard IHC Process




SCIZYS Erbium Kit



Preparation of the tissue section with SCIZYS Erbium Kit




Placement of t in SCIZYS slide

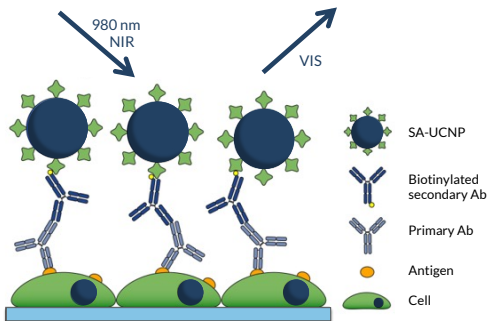


## Upconversion IHC-technology by Lumito

The UCNPs (Upconversion Nanoparticles) are illuminated by Near Infrared (NIR) light, emitting visible spectrum light. Lumito's NIR-imaging with UCNPs offers highly sensitive antibody detection and eliminates undesirable tissue autofluorescence.

## Invisible Labelling and Dual Modalities

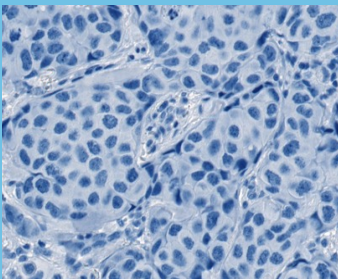

Brightfield  
Visible light

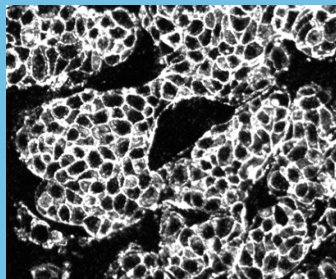


UCNP  
"Invisible" light

Images of dual emission modalities from separate optical channels.




the sample  
holder

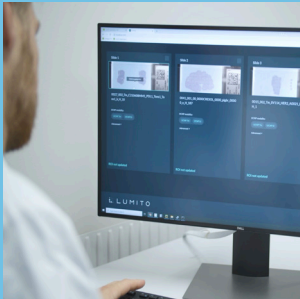


Automatic overview.  
Select your ROI area



Images delivered in high  
resolution for analysis in  
an image viewer

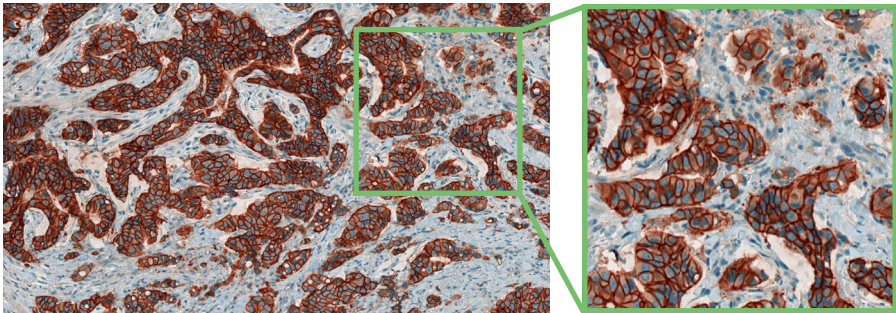

HER2 labelled with DAB

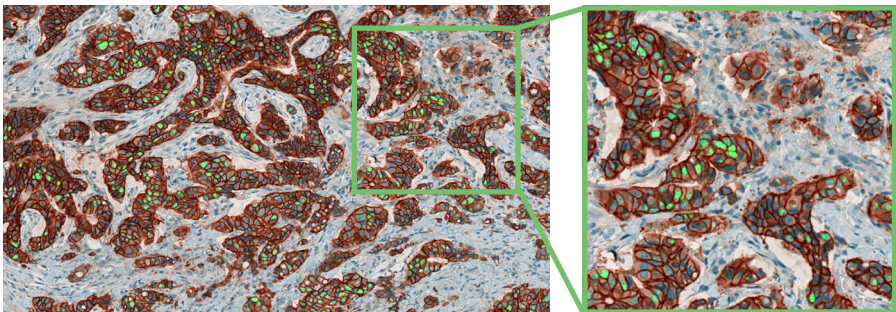

PR with UCNP as reporter

The images show that there is no interference with other chromogens and that the user extracts more information from each section by targeting different biomarkers with high accuracy and no overlap.



For Research Use Only

## **Elin Lindh**

Marketing Coordinator

Email: [info@lumito.se](mailto:info@lumito.se)

Tel: +46 (0)72 572 15 52



**Learn more!**

## **About Lumito**

Lumito is a spin-off from a research team at Lund University's Division of Atomic Physics and Laser Centre. Our mission is to offer a product that frees up resources, increases precision, and minimises human margins of error when diagnosing cancer indications.

Lumito AB, Mårtensstorget 5 SE-223 51 Lund Sweden  
[www.lumito.se/en/](http://www.lumito.se/en/)